


Single Use Hot Biopsy Forceps

FD-231C

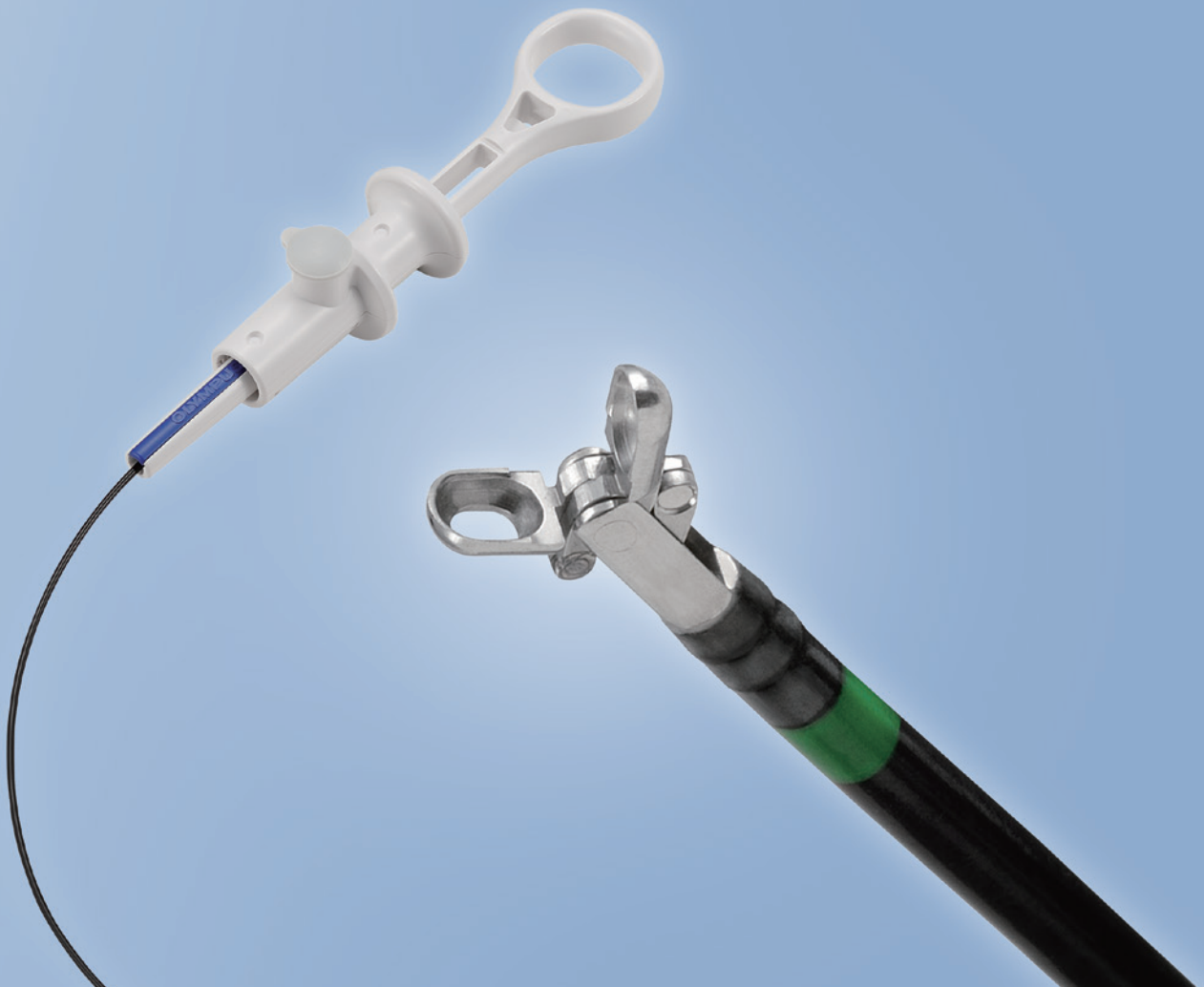


Single Use Hot Biopsy Forceps

Introducing the Hybrid Device-Disposable Hot Biopsy

Product Features


Hybrid device having both hot biopsy forceps and coagulation electrode function

- Used for electrosurgical collecting tissue, electrosurgical cauterization and hemostasis within the upper airways, tracheobronchial tree and urinary organs.
- Used for hemostasis same as coagulation electrode (CD-6C-1) when the cup of the forceps is closed.

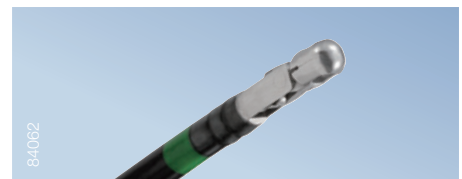

Marking for safety

Used as a reference in the endoscope field of view so that the metal portion and distal end of the endoscope are kept at a distance that allows safe activation of the electrosurgical unit.

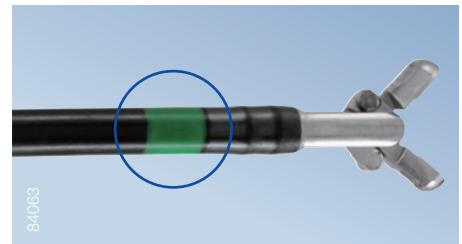


Compatible Olympus Electrosurgical Unit

ESG-150, ESG-300

Specifications

Model	Feature	Quantity	Working length	Sheath color	Minimum channel size
FD-231C	Oval fenestrated	1/box	1050 mm	Black	2.0 mm

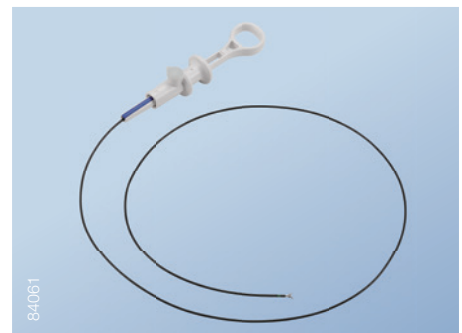

Olympus reserves the right of errors, modification and changes of the service and/or product offerings.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG
Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 14-18, 20097 Hamburg, Germany
Phone: +49 40 23773-0, Fax: +49 40 233765
www.olympus-europa.com